




CIDEX®

Activated Dialdehyde Solution

Excellent high-level disinfection from a 2.4% alkaline glutaraldehyde solution with a 14-day reuse life.

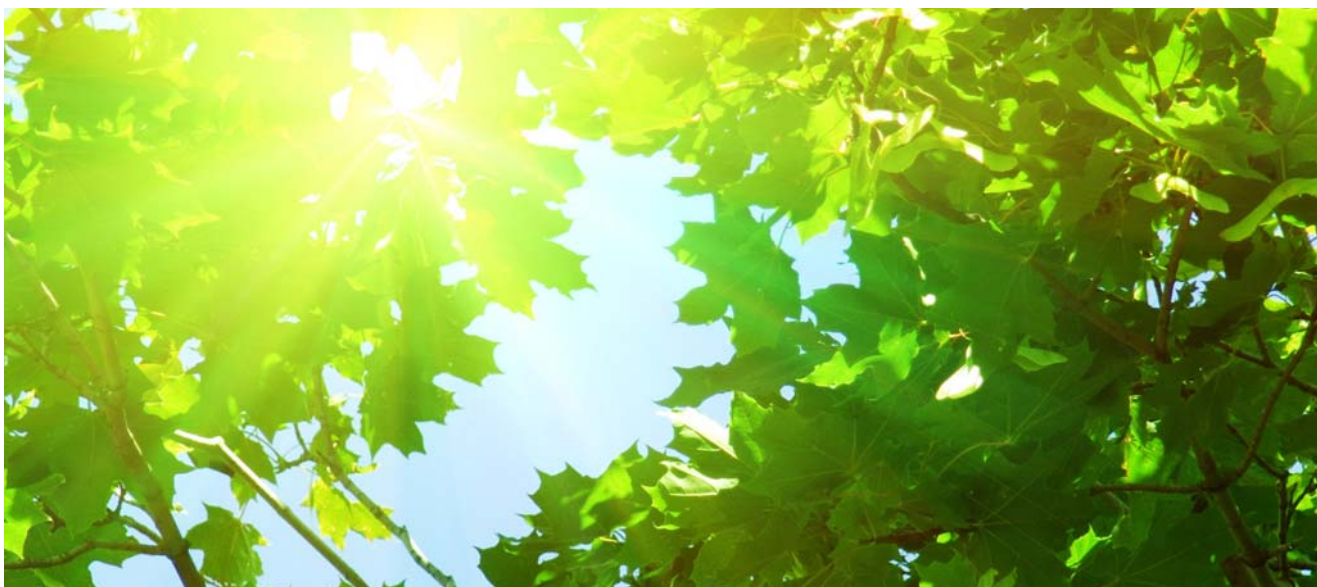

CHECK OUT THE NEW BOTTLE SIZE!

ASP is thinking green! The new CIDEX® Activated Dialdehyde Solution bottle uses less plastic, almost 20% more solution per bottle—and that means less waste, which is better for the environment. Plus, the IFU is available online at www.e-ifu.com, which means less paper and more trees! Product code 2266 now replaces product code 2250.

ADVANCED STERILIZATION PRODUCTS

Division of Ethicon, Inc.

a **Johnson & Johnson** company

Trusted by hospitals the world over, CIDEX® Activated Dialdehyde Solution provides rapid, cost-effective, high-level disinfection for a wide range of endoscopes and other medical instruments. CIDEX® Activated Dialdehyde Solution offers excellent tuberculocidal activity combined with the noncorrosive advantages of glutaraldehyde. CIDEX® Activated Dialdehyde Solution is a 2.4% alkaline glutaraldehyde solution that provides excellent microbicidal and anticorrosion properties vs. acid glutaraldehydes.¹

EXCELLENT COMPATIBILITY

- Solution is compatible with a wide variety of medical instruments.

SIMPLE ACTIVATION

- Activate the CIDEX® Activated Dialdehyde Solution by adding the entire contents of the Activator vial, which is included with the CIDEX® Solution container.

RELIABLE

- Accurate, easy-to-use test strips detect the minimum effective concentration (MEC).
- Activated solution offers up to 14-day reuse life.

ACHIEVES HIGH-LEVEL DISINFECTION

- Broad-spectrum antimicrobial activity destroys all vegetative forms of bacteria as well as 6 log reduction of *Mycobacterium tuberculosis* in 45 minutes at 25° C.
- Effective in the presence of 2% organic soil.

SIMPLIFIES DISINFECTION

- After activating, just pour activated solution into compatible tray or automatic endoscope reprocessor and completely immerse instruments.
- Conveniently discarded down hospital and office drains in accordance with federal, state and local regulations.

ORDERING INFORMATION

REORDER#	DESCRIPTION	CASE CONTENTS
2266	Container of CIDEX® Activated Dialdehyde Solution	1.24 gallons (4.7L), 4 bottles
2920	CIDEX® Activated Dialdehyde Solution Test Strips	60 Strips/bottle, 2 Bottles

Please read and follow the CIDEX® Activated Dialdehyde Solution instructions for use prior to using this product for complete usage information, including contraindications and warnings.



1. Rutala W. APIC guidelines for selection and use of disinfectants. Am J Inf Cont. 1996;24:313-342.
© Ethicon, Inc. 2011. All rights reserved.